



PD Care Laser

FOR PATIENTS



About SYMBYX & Light (PBM) Therapy

What We Do?

SYMBYX is an Australian medical technology company developing light therapies to manage previously intractable, chronic diseases and painful inflammatory conditions. A key focus for SYMBYX is using light therapy (photobiomodulation or PBM) to reduce Parkinson's disease symptoms. We treat primarily via the gut microbiome because there is now a well-known and researched linkage between the gut and the brain (known as the Gut-Brain-Axis).

Light Therapy and How it Works?

PBM is not a new area of medicine. In fact, the 1903 Nobel Prize in medicine was won by a Danish doctor for using ultraviolet blue and red light to treat smallpox and TB scars. Today, there are over 4,000 peer-reviewed scientific studies reporting on light's safe, protective and restorative properties for a range of health problems.

Applying an optimal dosage of light (which considers wavelength, application time and intensity or strength) is critical – application above or below this dose will have diminished effect. SYMBYX lasers deliver a specifically calibrated light energy doses relying on proprietary, super-pulsing technology. This allows for the delivery of optimal dosages without thermal (heating) effect. Our lasers are also non-invasive, painless to use and portable.



Unlocking the
healing power of
light

SYMBYX light therapy can relieve Parkinson's symptoms



Why is Light Relevant to Parkinson's disease?

In Parkinson's disease it is not only the loss of dopamine producing brain cells but also mitochondrial dysfunction (a part of every living cell which produces the energy required for daily life) that causes both movement symptoms and other non-motor symptoms — such as diminished balance, cognitive ability, sense of smell, problems with sleeping, depression, and loss of motivation. Light therapy energises cell mitochondria and can help with the loss of cellular energy common in Parkinson's disease.

Most people don't realise that their colon has the capacity to produce both dopamine and serotonin – critical neurotransmitters required for healthy brain function. PBM to the gut reduces inflammation, promotes the metabolism of short-chain fatty acids (SCFAs), which in turn influences the production of more dopamine and serotonin. This process delivers therapeutic benefits to people with Parkinson's.

People with Parkinson's are reporting improvements in mobility, cognition, dynamic balance, sleep, motivation and fine motor skills after using the Symbyx PDCare laser. BMC Neurology (July 2021) noted "no adverse side effects" and that PBM was shown to be a "safe and potentially effective treatment for Parkinson's disease symptoms. Improvements were maintained for as long as treatment continued, for up to one year in a neurodegenerative disease where decline is typically expected".

Our painless, easy-to-use SYMBYX lasers are CE marked and ARTG listed.



PDCCare 904 nm Laser

All SYMBYX devices are approved for home usage. SYMBYX lasers are high-quality, European manufactured for superior precision, increased reliability, and consistency over time. It is recommended that the diodes are re-calibrated every 3-5 years for optimal performance. Lasers significantly reduce treatment times versus simple light emitting diodes (LEDs).



Each laser is delivered in a package that includes:

- PDCCare 904 nm Laser by SYMBYX
- 6 x AAA batteries
- Instruction manual
- Protective sleeve for device



1. Laser aperture/Laser diodes

The laser emitting diodes are placed in a slight arch to improve contact between the laser apertures and the body in order to optimise penetration of light into the tissue. There are two diodes.

2. Powder painted aluminum front

The powder coated aluminum front and the anti-reflective coated lenses ensure long-term use and allow cleaning without deterioration.

3. Feedback/ LED light indicator

A semi-transparent ring emits a green light indicating that laser emission is activated. Feedback is also given by visible and audible signals when sufficient laser energy has been delivered and when it is time to change the treatment point. This feedback occurs every 60 seconds.

Laser emission cannot occur accidentally. The operator must first turn on the device into stand-by mode and also then lightly but firmly press the device directly onto clean, dry skin to commence transmission.

4. Concave shape for better grip

The laser surface has a concave shape along the length of the device, which provides a better grip, is more comfortable to hold and reduces the risk of dropping the device.

5. Screw for replacing batteries

The laser is internally powered by 6 AAA-batteries and has a battery-life of approximately 3-4 treatment hours. Remember to switch it completely off between uses. Rechargeable batteries may be used.

6. Power on/off

Accessible button on the back to turn the laser on and off out of stand-by mode.

Treatment Protocol

The PDCARE Laser is to be used to treat the neck and abdomen.

Time to complete: 20 minutes

18 minutes on stomach 2 minutes on neck at base of the skull.

Frequency: 3x per week

Turn the laser on with the button on the base (6) of the device - you will hear 4 short beeps and see a flashing light (3).


Press it onto the spot (2) to be treated with enough pressure to turn it on. You will hear one beep.

Thereafter, the laser will beep every 1 minute. You should hold it on the 9 spots marked on the abdomen for 2 full beeps (2 minutes x 9 spots = 18 minutes total). **Treat over your bare abdomen** in a rough grid pattern as shown in the diagram. **The laser needs to make direct contact with the skin.** The top 3 spots are just under the ribs and the bottom 3 spots are just above the pelvis.

After completing the stomach routine, press the laser against the neck at the base of the skull for 2 beeps (2 minutes). Only do this on days you use the laser. Each complete laser session is therefore 20 minutes long. Turn off the device using the button on the base (6) of the laser.

We encourage you to seek a follow up appointment with your Clinician after a 30-day trial for a review before making your decision to keep or return the laser.

Also, if you any questions during the trial period contact your clinician or email SYMBYX on clinicalsupport@symbyxbiome.com

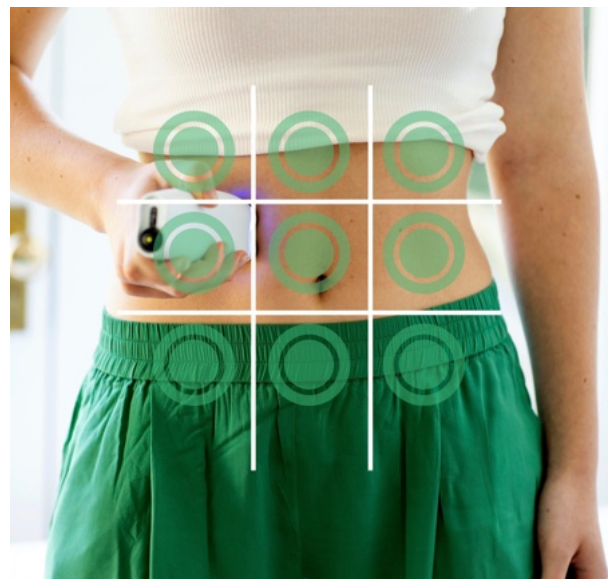

PDCare Laser: Clinician Supported Trial

Symptoms of Parkinson's disease vary from person to person and because of this, treatments that work for one person may not work the same way for another. Understandably, this means deciding to invest in a new treatment therapy can be a big decision.

To help you make this decision, you can trial the SYMBYX PDCare laser under the care of your Clinician in their clinic before deciding whether to buy your own. Your Clinician will show you how to use the laser and will conduct some assessments. These assessments, along with the advice of your Clinician, will help you decide whether continuing with the treatment at home with your own laser is the best course of action for you.

Our clinical research shows that the the best outcomes for people with Parkinson's are achieved when the PDCare laser is used three times a week in accordance with the SYMBYX protocol. For this reason, the therapy is ideally self- administered in the comfort of your own home.

You can order your SYMBYX laser online by following these steps:

<https://symbyxbiome.com/> > Our products > add PDCare 904nm Laser to cart > add discount code (provided by your clinician) > fill out shipping address > add billing details > finalise purchase



Enter your shipping details (device will be sent to you) > finalise purchase.

If you decide to purchase, use the exclusive code your clinician gives you so you qualify for a 10% discount off our RRP.

Please reach out to us to learn more

Website: <https://symbyxbiome.com>

General enquiries: info@symbyxbiome.com

Clinical questions: clinicalsupport@symbyxbiome.com

